
SALVAFAST

TRANSFER ***RESTRAINT SYSTEM***

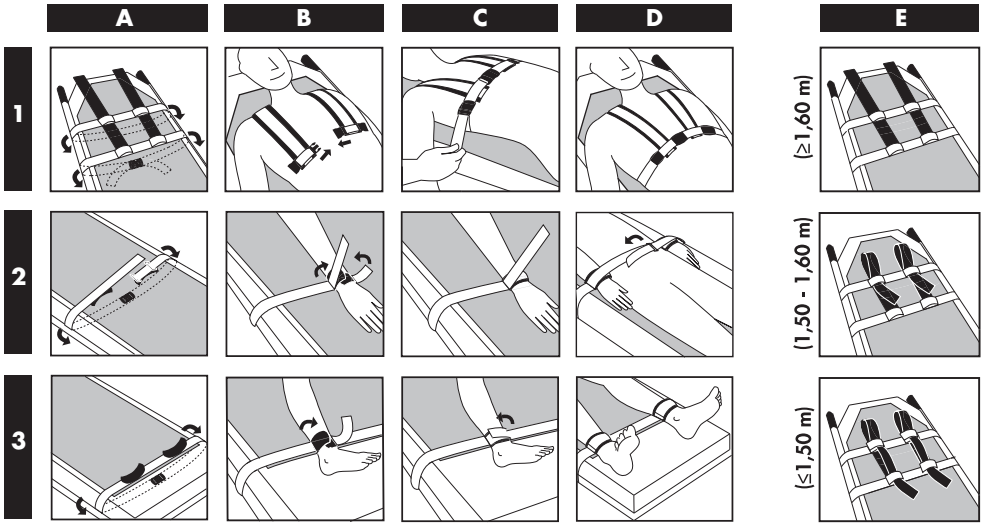MEDI CARE SYSTEM, S.L.U.
Abat Escarré, 15
08302 Mataró - SPAIN
info@medisystemsl.com
www.medicaresystem.es

CE MD

MEDICARE SYSTEM

ESPAÑOL

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ENGLISH

APPLICATION INSTRUCTIONS

• Unfold the product over the stretcher and place the round side of the mat on the top.
• Tighten the blue straps under the stretcher. **A.**
• Lay down the patient over the mat and fasten the ankle bands. **3.**
• Place the harness over the patient's shoulders. Use the loop straps to adjust the harness to the patient's height. **E**
• Check that there is no looseness between the shoulder and harness and fasten the buckles over the patient's chest.**1.**
• Fasten both wristbands and insert the remaining band through the adjustable rings. Tighten the strap over the patient's waist and fasten the strap over itself using the Velcro® fasteners. **2.**
• The handles are identified in orange colour under the blue straps. For patient transfer, unfasten the blue straps buckles and use the handles to transfer the patient from a stretcher to a bed.

SAFETY

1. Follow the instructions of the restraint protocol of your health centre and the laws of your country.
2. Fasten all the items specified above, harness, wrist and ankle bands.**1-2-3.**
3. Salvafast Transfer is a restraint system for patient's transport. Patients must be monitored all time while being restrained.

WARNING

• DO NOT USE the product if you do not comply with ALL these points.
• Improper use of the product may cause serious injury.

MAINTENANCE

• Clean the product after use with a damp cloth and mild soap. In case of disinfection, follow the disinfection guidelines of your health-care facility. Use a surface virucidal disinfectant in spray or wipes. Aqueous solutions of ethanol or isopropanol. The product can be washed at 60°C if required using specific disinfectant detergent. However, it is recommended to hand wash the product for a longer life cycle.
• Salvafast products are not fireproof.
• Fasten Velcro® fasteners and buckles before machine washing.

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